



# HOVERMATT®

## Single Patient Use

The number one choice for lateral patient transfers now offers the added convenience of disposability in the Single-Patient Use HoverMatt® transfer mattress. The inflatable transfer mattress works just like the reusable HoverMatt® air transfer system, using a cushion of air that enables caregivers to safely transfer patients without lifting or straining. By eliminating injuries related to lateral transfers and repositioning, the Single-Patient Use HoverMatt® can help improve staff retention while meeting legislative guidelines for safe patient handling.




Variable Speed  
HoverMatt® Air Supply

### HoverMatt® Single Patient Use Benefits:

- Improves patient and caregiver safety and comfort during transfers
- Reduces skin shear and bruising
- Eliminates reprocessing (no laundering)
- Ideal for use in any patient care department
- Radiolucent/MRI compatible
- Available in multiple styles and widths to accommodate all patient sizes
- Significantly reduces the possibility of cross contamination
- Reduces costs and logistics associated with reprocessing
- Variable speed air supply provides rapid inflation
- Weight Limit 544kg

**HoverMatt Single Patient Use**

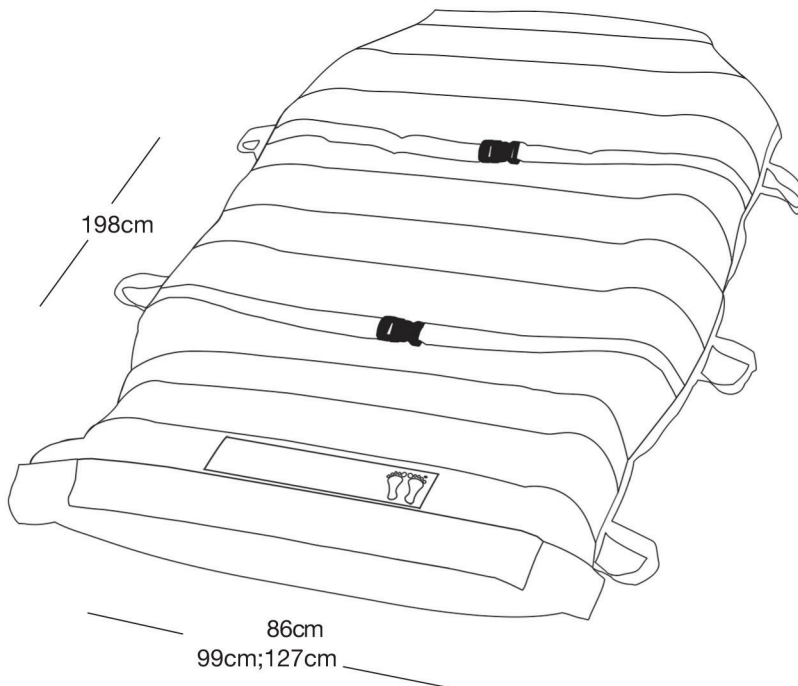

## Product Specifications

### Material Specs

Sewn Construction

Top: Non-Woven Polypropylene Fibre

Bottom: Nylon Fabric

Latex free

Water Soluble Perimeter thread

(120°) - Thread in sewn handles is NOT Water Soluble

### Available in 5 Sizes

86cm W x 198cm L (MHHT-HMSPU-34)

99cm W x 198cm L (MHHT-HMSPU-39)

127cm W x 198cm L (MHHT-HMSPU-50)

86cm W x 115cm L (MHHT-HMSPU-34MI)

99cm W x 115cm L (MHHT-HMSPU-39MI)

### MHHT-HMSPU-34, MHHT-HMSPU-39, MHHT-HMSPU-34MI, MHHT-HMSPU-39MI

10 per box

### MHHT-HMSPU-50

5 per box

### Compatible with HoverTech Air Supplies

MHHT-AIR2300 (Variable Speed)

*If a HoverMatt® SPU is laundered, it will no longer inflate when connected to the air supply. The "DO NOT WASH" label at the foot-end of the product will dissolve and state "DISCARD IMMEDIATELY".*

For more information, contact your local HoverTech distributor, or Statina Healthcare International on 1300 791 404.

### Distributed By: